# S · E · Y · S · S · O

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Before using this product, please read carefully this user's manual and keep it properly for later reference.

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# English

<b>1</b>	Attention	4
<b>2</b>	Articles in the packaging box	6
<b>3</b>	How to charge your SEYSSO sonic toothbrush	7
<b>4</b>	To carry the SEYSSO sonic toothbrush with you	8
<b>5</b>	How to use your SEYSSO sonic toothbrush	9
<b>6</b>	How to clean your SEYSSO sonic toothbrush	11
<b>7</b>	How to store your SEYSSO sonic toothbrush	12
<b>8</b>	Environmental protection	12
<b>9</b>	To replace brush heads	13
<b>10</b>	Warranty and services	13

# 1 Attention

**Before using this product, please read carefully this user's manual and keep it properly for later reference.**

## **Danger:**

- Keep the inserted charging base away from water. Do not place it above or nearby any container with water (e.g. bathtubs or washbasins). Do not place the charging base or rechargeable travel storage box in water or any other liquid. After cleaning the product and before connecting it to the power, please ensure that the charging base or the rechargeable travel storage box is thoroughly dry.

## **Warning:**

- Before connecting the product to the power, please check whether the voltage marked on the product agrees with the local supply voltage.
- If the charging base or the rechargeable travel storage box has any damage, please stop using the product. If the charging base, rechargeable travel storage box, power line or plug is damaged, please be sure to use a new one of the original model to replace the damaged one to avoid danger.
- If the product has any damage (e.g. brush head or toothbrush handle), please stop using it. The product contains some parts that cannot be replaced by user. If the product is damaged, please contact the customer service center in your country/region (please refer to the chapter "Warranty and Services").
- Please do not use charging base or rechargeable travel storage box in the open air or nearby an article with high-temperature surface.
- This product is not applicable to those who are disabled in body, sense or spirit or lack relevant experience and knowledge (including children) unless their guardians supervise or instruct them how to use this product.
- Please look after your kids and ensure that they do not play with this product.


## Attention:

- If you had oral or gum operation in the past two months, please consult your doctor before using this product.
- If your gum bleeds excessively after using this product or still bleeds after using it one week later, please consult your doctor.
- SEYSSO sonic toothbrush accords with the safety, standard for electromagnetic equipment. If you have any question about cardiac pacemakers or other implanted devices, please consult your doctor or the manufacturers of the implanted devices before using SEYSSO sonic toothbrush.
- This product is designed only for cleaning teeth, gum and tongue. Please do not use it for other purposes. If you feel not comfortable or painful, please stop using this product and contact your doctor.
- If you have any medical problems, please consult your doctor before using this SEYSSO sonic toothbrush.
- Please stop using deformed brush heads or those with curved bristles. Please replace your brush head every three months. If the brush head is worn, please replace it ahead of time.
- Please do not clean your brush head in a dish-washing machine.
- Please do not use your brush heads not recommended by the manufacturer.
- Please give guidance to children for tooth brushing until they form good tooth brushing habits.
- If the toothpaste you use contains peroxide, sodium bicarbonate or other bicarbonate (commonly seen ingredients in whitening toothpaste), it is suggested that you clean your brush head thoroughly with soap-suds after using it every time.

## EMF:

This product accords with all relevant standards for EMF. According to current scientific proof, if you use this product correctly in accordance with this user's manual, this product is safe.

## 2 Articles in the packaging box

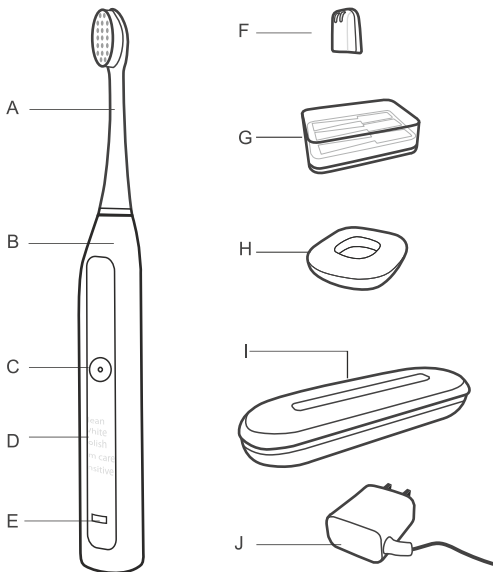- A. Brush head
- B. Toothbrush holder
- C. Power switch button
- D. Tooth brushing mode
- E. Charging indicator

- F. Travel brush head protective cap
- G. Brush head storage box
- H. Inductive charging base
- I. Rechargeable travel storage box
- J. MicroUSB charger


### 3 How to charge your SEYSSO sonic toothbrush

Insert the toothbrush handle into the charging base to charge it. Electricity quality is showed on the LED electricity quality display. The fast blinking red indicator means low battery (for less than three times). The charging indicator should blink through the charging process and turn to continuous light when battery is fully charged. It will take about 15 hours to fully charge the handle, and it can be used up to 70 days straight.

Result depends on a chosen brushing mode and regards using twice a day:


- a. **Clean:**>57 days;
- b. **White:**> 25 days;
- c. **Polish:**> 45 days;
- d. **Gum Care:**> 58 days;
- e. **Sensitive:**> 70 days.

If the battery is low, the electricity quality display LED light will flash and the speed of motor will reduce. When the battery runs out, the motor will stop working. It needs at least 30 minutes to change for reboot and hour later it can start to work. For convenient use, it is suggested that you use this product after it is fully charged.
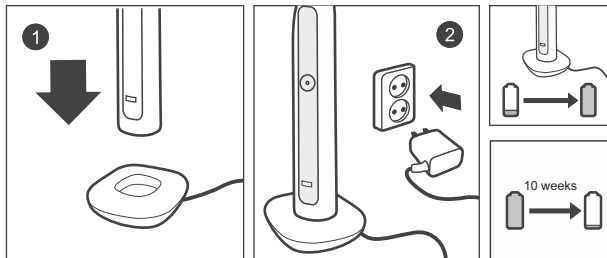
## 4 To carry your SEYSSO sonic toothbrush with you

Put the toothbrush into rechargeable travel storage box. Connect the microUSB interface with the rechargeable travel storage box. The charging indicator should blink through the charging process and turn to continuous light when battery is fully charged.

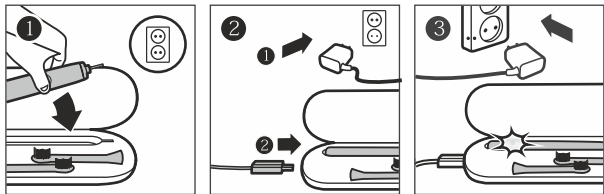
## 5 How to use your SEYSSO sonic toothbrush

**Tooth brushing:** gently and slowly move the brush head between your teeth, so that the long brush bristles can clean your gum lines. Please do not brush your teeth with might like using a manual toothbrush. To give full play to the function of SEYSSO sonic toothbrush, you should exert a bit force to make SEYSSO sonic toothbrush serve you.


**Intelligent zoning:** it divides your oral cavity into four areas and it will remind you to change to another area per 30 seconds

**Two-minute timing function:** when two-minute scientific tooth brushing time is reached, the electric toothbrush will automatically stop working.

**Intelligent memory chip:** it will remember your using habit and will switch to the gear you used last time when it is started again.

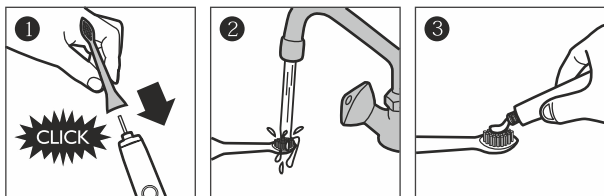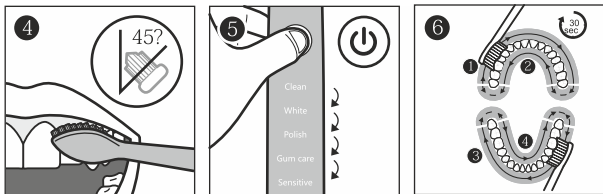Five modes of SEYSSO sonic toothbrush:


- Clean** 82.000 mpm\* – Efficiently removing dental plaque (*2min*)
- White** 96.000 mpm – Getting rid of surface stains and whitening teeth (*2min*)
- Polish** 90.000 mpm – Brightening and polishing front teeth (*1min*)
- Gum Care** 70.000 mpm – Gently stimulating and massaging gum and improving gum health (*3min*)
- Sensitive** 62.000 mpm – An extremely gentle mode for sensitive teeth and gums (*2min*)

\*Mpm: movements per minute

## 6 To clean your SEYSSO sonic toothbrush


Please regularly clean your SEYSSO sonic toothbrush, wash the toothbrush handle and head. If necessary, please use soft cloth and mild soap to clean its surface. Please do not use sharp articles to clean it. After each use remove the brush head. Carefully clean the handle top and the pin to remove the toothpaste sediment.

Please do not clean the brush head, handle, charging base or rechargeable travel storage box in dish-washing machine.
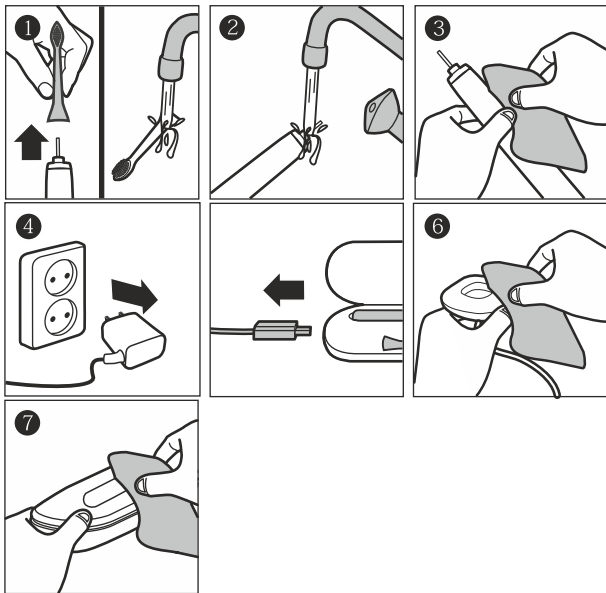
## 7 How to store your SEYSSO sonic toothbrush

If you will not use your SEYSSO sonic toothbrush for a long time, please disconnect the power of the charging base, clean all parts and store the toothbrush in a cool and dry place out of sunshine.
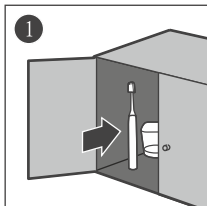
## 8 Environmental protection

If you want to discard this product, please do not place it together with domestic waste but send it to the officially designated recycling center and dispose the battery at the officially designed battery recycling center. Doing so can benefit environmental protection.

## 9 Replace brush head

SEYSSO warranty will be invalid if it is found that the handle is damaged because non-SEYSSO replacement brush heads are used.

It is not recommended to use non-SEYSSO replacement brush heads, because SEYSSO cannot control the quality of other replacement brush heads.

SEYSSO cannot guarantee that non-SEYSSO replacement brush heads fit the handle.

All SEYSSO replacement brush heads are printed with the logo of SEYSSO and meet the high-quality standard of SEYSSO.

## 10 Warranty and services

You will enjoy two-year quality warranty for this product since the day of purchasing. In the warranty period, we will be responsible for repairing or replacing the complete device if malfunction is caused by any design defect, but the warranty does not cover the damages caused by inappropriate use of the product.

### Quality warranty – restrictive clause

Quality warranty is only applicable to the products manufactured by SEYSSO.

Official website: [www.seysson.com](http://www.seysson.com)

## Explanation

Please carefully fill in the warranty card and keep it appropriately for later service. You can enjoy two-year free warranty services for the product since the day of purchasing. You cannot enjoy free services if any of the following circumstances happens:

1. The product is damaged because of inappropriate use or improper storage;
2. The product has been disassembled or repaired by others rather than SEYSSO's special repair centers;
3. You do not have effective proof of purchase;
4. Damages are caused by instable voltage or voltage beyond the specialized range;
5. The product is not in the warranty period or beyond the coverage of warranty.

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