S · E · Y · S · S · O

Edge

INSTRUKCJA DOTYCZY MODELI O NUMERACH
/ THE INSTRUCTION APPLIES TO MODELS
WITH THE FOLLOWING NUMBERS / TENTO
NÁVOD PLATÍ PRO MODELY S NÁSLEDUJÍCIMI
ČÍSLY / TENTO NÁVOD PLATÍ PRE MODELY S
NÁSLEDUJÚCIMI ČÍSLAMI:
SE.WF.010.BLK.1, SE.WF.010.WHT.1

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
NOZZLE DESCRIPTIONS:

- Standard nozzle (x2): For daily cleaning.
- Tongue cleaner nozzle (x1): Removes bacteria and ensures fresh breath.
- Plaque Seek nozzle (x1): Designed for cleaning around implants, crowns, and bridges.
- Orthodontic nozzle (x1): Removes food debris from braces.
- Pik Pocket nozzle (x1): Flushes and massages gum pockets.

PRODUCT DESCRIPTION

1. Irrigator nozzle
2. Rotating nozzle holder
3. OLED display
4. Power button
5. Pressure adjustment button "+"
6. Mode selection button
7. Pressure adjustment button "-"
8. Water tank
9. Nozzle release button
10. Drain hole
11. USB charging port
12. Temperature sensor
13. Water supply tube
14. Gravity ball valve
15. UVC disinfection lamp

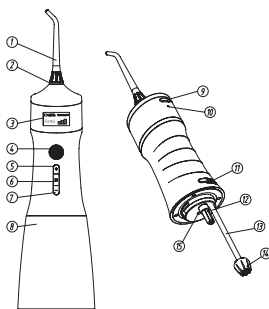Note: During the first use, there may be a small amount of distilled water in the tank, which was used for production testing. It is completely safe for health, and there is no cause for concern.

WARNINGS! READ BEFORE USE! KEEP THIS MANUAL..

When using electrical devices, especially in the presence of children, always follow basic safety precautions, including the guidelines below.

- For household use only. 
- The irrigator has a waterproof casing. It is safe and designed for use in the bathroom. The device has an IPX7 protection rating, meaning it is resistant to accidental immersion in water up to 1 meter deep for up to 30 minutes.
- Important: The irrigator is not designed for underwater operation!
- The product is not intended for self-repair.
- Keep the charging cable in a dry and safe place.
- The device should be charged only with the cable included in the set.
- Straighten the cable before use.
- Do not place the device in fire, hot or corrosive liquids, or drop it into them.
- Do not reach for the device if it has fallen into water. Immediately disconnect it from power.
- In case of malfunction, do not disassemble the device or replace any components, as there is a risk of electric shock. Contact SEYSSO service.
- If the power cord is damaged, it must be replaced with a new one. Contact SEYSSO service.
- If dropped, replace the nozzle before the next use, even if no visible damage is present.
- Do not disassemble the device except when disposing of it. The battery must be removed before disposal.
- Do not dismantle the device or replace its components unless you are a specialist.
- Usage Precautions:
 - Do not fill the water tank with water hotter than 40°C (104°F).
 - Do not turn on the device without an attached nozzle.
 - Keep the irrigator nozzle out of children's reach—risk of choking.
 - Mouthwash liquids may damage the device. Do not pour mouthwash into the water tank.
 - If you experience pain, discomfort, or concerns about the device's impact on your health, consult a doctor.
 - If you have had oral surgery within the last two months, consult a doctor before use.
 - Do not press the nozzle too hard against gums or teeth.
 - The device is designed for cleaning interdental spaces, gums, and tongue.
 - Do not direct the water stream under the tongue, into the ear, nose, or

other sensitive areas. The device generates pressure that may cause serious injury. Use only as intended.

- Use SEYSSO Edge Pocket nozzles only on the Soft mode. High pressure may damage the nozzle.
- Do not connect the device to a power system with a voltage different from that specified on the device or charger.
- Do not drop or insert foreign objects into the device's openings or tubing.

DANGER WARNING

To reduce the risk of electric shock:

- Always unplug the device from the power source after use.
- Do not use while bathing or showering.
- For detailed instructions on storage, usage, and included accessories, contact SEYSSO.
- To reduce the risk of fire, short circuit, or injury:
- Unplug the charger from the power source after charging is complete.
- This product is not a toy. Exercise caution when the irrigator is used by or near children.
- Use only as intended. Only use accessories recommended by the manufacturer.
- Do not use the charger if the power cable or plug is damaged. If the irrigator malfunctions after being dropped, left outdoors, or submerged in water, return it to the place of purchase.
- If the power cable (and/or plug) falls into water, do not attempt to retrieve it. Immediately unplug it from the power outlet.
- Do not pull, pinch, or wrap the power cord around sharp edges. Avoid closing it in doors or exposing it to heat.
- Do not turn on the irrigator while charging. Keep all components away from heated surfaces.
- Do not unplug the charger by pulling the cable. Always grasp the plug directly.
- Do not handle the plug with wet hands.
- Never use the product if you have trouble concentrating or feel drowsy.
- Do not use the device outdoors or in areas exposed to aerosols or pure oxygen.
- Plug the device into a suitable power outlet and ensure compatibility with its specifications.
- Children & Special Considerations:
- Not recommended for children under 8 years old.
- Children and individuals with sensory, motor, or cognitive impairments should only use the device under supervision or after receiving proper training on safe usage.

ADDITIONAL NOTES

- Do not share your nozzle with other users.
- Do not tightly coil the power cable.
- Unplug the charging cable from the outlet if it will not be used for an extended period.
- Do not use the device if you have an open wound on your tongue or in your mouth.
- Do not use iodine or water-insoluble concentrated essential oils in this device. These substances may reduce performance, shorten product lifespan, and void the warranty.
- Children over 8 years old should be supervised to ensure they do not play with the device.
- Never leave the device connected to power unattended.
- If the product is opened or disassembled for ANY reason, the warranty is void.
- To prevent damage, do not place the device or its parts in the dishwasher, unless otherwise specified by the manufacturer.
- Use the device only as instructed or as recommended by your dentist.
- Use only SEYSSO-approved nozzles and accessories.
- Fill the water tank with water only.
- Remove all oral jewelry before use.
- If your doctor or cardiologist has recommended taking antibiotics before dental procedures, consult your dentist before using the device.
- Contact your dentist if you experience excessive bleeding after use or if bleeding persists for more than a week.
- Stop using the device and consult a doctor or dentist if you experience discomfort or pain.
- The device meets electromagnetic safety standards. If you have a pacemaker or other implanted medical device, consult your doctor or the device manufacturer before using it.

GETTING STARTED

Powering the Device

1. Fully charge the device before first use.
2. Use the original USB cable included in the set.
3. Connect the USB end to a standard DC 5V 1A port and the Type-C end to the device's charging port.

Charging the Device


1. During charging, a charging indicator will appear on the display (see Image 2).
2. After 4-5 hours, when fully charged, the screen will display „100% Com-

pleted" (see Image 3).

3. When the „Low Power" message appears, recharge the device (see Image 4).

4. While charging:

- Keep the charging port dry and do not rinse the device with water.
- Place the device upright on a desk; do not lay it flat.
- Do not cover the device with a towel or wrap it while charging.

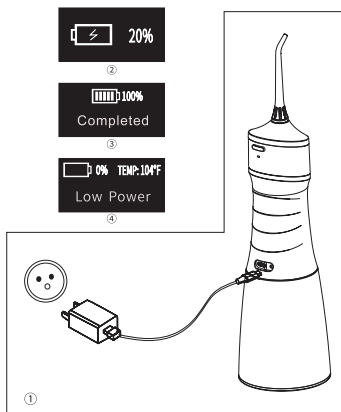Inserting and Removing Nozzles

1. Inserting the Nozzle:


- Insert the nozzle into the central opening until you hear a click (see Image 5).

2. Removing the Nozzle:

- Press the nozzle release button (see Image 6).

3. Important Notice!

- Do not press the nozzle release button while the device is operating, as the nozzle may eject forcefully, potentially causing injury.

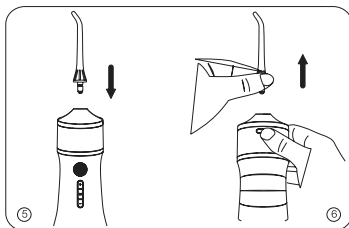Operating the Device



1. Turning on the Device

- Press and hold the ON/OFF button for 2 seconds.
- The device will turn on and enter standby mode.
- The screen will display battery status, water temperature, mode, and pressure level (see Image 7).

2. Changing the Cleaning Mode

- Short press the „M“ button to switch between cleaning modes (see Image 8).
- Mode cycle:
- Strong → Normal → Soft → Pulse

⑧

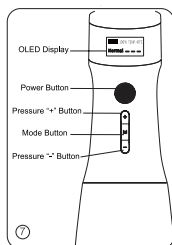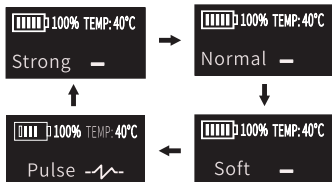
3. Changing the Temperature Unit

- While in standby mode, press the „+“ and „-“ buttons simultaneously.
- The device will switch the temperature unit between Celsius („°C“) and Fahrenheit („°F“).


4. Cleaning Mode Descriptions

- Strong Mode: Powerful water stream, ideal for deep cleaning and removing stubborn deposits.
- Normal Mode: Medium power, suitable for daily use for most users.
- Soft Mode: Gentle water stream, recommended for sensitive teeth, new users, and children over 8 years old.
- Pulse Mode: Alternating cycle of strong and weak streams, gently massages the gums, recommended for mild gum inflammation.

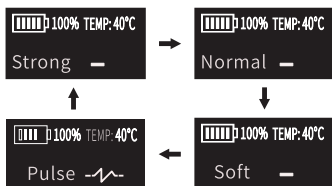5. Turning the Device On and Off

- After selecting a mode, short press the ON/OFF button to start the device.
- Short press the ON/OFF button again to pause the device.

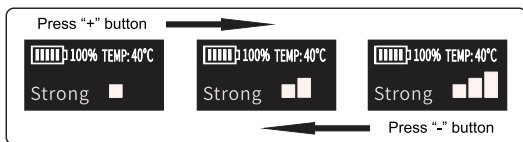6. Adjusting the Pressure Level

- Strong, Normal, and Soft modes have three adjustable pressure levels.
- Example: In Strong mode, pressing the "+" or "-" button will increase or decrease the pressure level.
- The display shows the selected pressure level (see illustration).

7. UVC Sterilization Mode

- When the device is turned off, press and hold the „M“ button for a few seconds to activate UVC mode.
- Short press the ON/OFF button, and the UVC lamp will start a 120-second sterilization process.
- During sterilization, the screen will enter standby mode.

8. Completing the Sterilization Process

- After 120 seconds, the sterilization process will end.
- The display will return to standby mode automatically after 10 seconds.
- Alternatively, short pressing the ON/OFF button will immediately restore standby mode.


9. Turning Off the Device

- Press and hold the ON/OFF button for 0.5 seconds to turn off the device.
- If the device remains in standby mode for 35 seconds without activity, it will turn off automatically.

10. Accurate Temperature Measurement

- To obtain precise temperature readings, ensure that the water level in the tank is above the sensor marker, so the temperature sensor is fully submerged (see illustration 9).

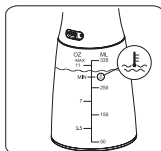11. Starting Oral Cleaning

Removing the Water Tank

1. Hold the device body and the water tank.
2. Twist the tank counterclockwise, then pull it downward to remove it (see illustration 10).

Filling the Water Tank

3. Fill the tank with lukewarm or filtered water (see illustration 11).

Reattaching the Water Tank

4. To reattach the tank, twist it clockwise until it clicks into place.

Turning on the Device

5. Select the desired mode and adjust the pressure level.



Using the Device

6. Insert the nozzle and position it at a 90° angle to the gum line or between the teeth.

7. Short press the ON/OFF button to start the water stream.

8. Move the nozzle along the gum line, gently shifting it back and forth.

9. Short press the ON/OFF button again to pause the device (see illustrations 12 and 13).

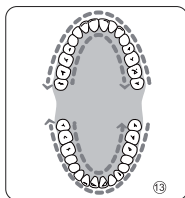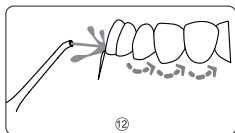12. Usage Tips



- If you wear braces or dentures, pay extra attention to cleaning interdental spaces.
- If you have sensitive teeth, start with the Soft mode and gradually adjust the pressure level according to your comfort.
- Do not tilt the device excessively during use, as this may cause the gravity ball valve to be out of the water, preventing proper water suction.

After Use & Maintenance

I. Cleaning the Water Tank

1. After removing the tank, use a bottle brush to remove any deposits inside the tank (see illustration 14).

2. Rinse the water tube and gravity ball valve or sterilize them using warm water (see illustration 15).

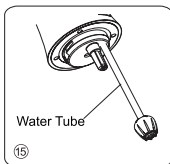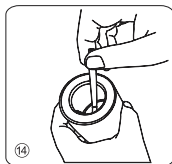II. Device Maintenance

1. Turn off the device after each use.
2. If the device will not be used for an extended period:
 - Empty the water tank completely.
 - Wash the tank and allow it to air dry before storing.
3. For hygiene reasons, replace the nozzle every 3 to 6 months, depending on the type of nozzle used.
4. Clean the device with clean water or warm water below 35°C (95°F).
5. Regularly clean the filter mesh at the end of the gravity ball valve to remove any debris.

Device Usage & Environmental Protection

- Do not dispose of the SEYSSO Edge Nova irrigator with household waste.
- Recycle it at an electronic waste disposal and recycling facility.

PRODUCT SPECIFICATION

Specification	Szczegóły
Model	SE.WF.010.BLK.1/SE.WF.010.WHT.1
Battery	1800 mAh Lithium-ion
Input	DC 5V = 1A
Pulse Frequency	2000-2800 pulses/min
Water Pressure	2,75-10,3 Bar
Tank Capacity	335 ml / 11.33 oz
Waterproof Rating	IPX7
Dimensions	12,2 x 3,2 x 3,5 inches

Important: Using the Correct Power Adapter

Please ensure that the power adapter used matches the output parameters specified on your irrigator's label. Using an adapter with incompatible specifications is strictly prohibited, as it may cause malfunctions and damage the device. Failure to follow this rule may void the warranty. Take care of your device by only using recommended accessories and power adapters.

TROUBLESHOOTING TABLE

Problem	Cause	Solution
Water leakage between the nozzle and base	Nozzle is not properly inserted	1. Remove the nozzle and reinsert it. 2. Replace it with a new nozzle.
No water stream	Insufficient water in the tank	Refill the water tank.
	Water tube is bent or gravity ball valve is clogged	Detach the water tank and check the tube and gravity ball valve.
	Air is trapped in the pump	Remove the nozzle and put a few drops of water into the nozzle opening.
Device does not turn on	Low battery level	Charge the device.
	ON/OFF button was not held long enough	Follow the user manual instructions.
	Internal electrical failure	Contact customer service for repair.
Screen does not work	Loose connection in the charging port	Replace the charging cable or contact customer service.
	Damaged charging adapter	Replace the adapter with a new one.
	Internal electrical failure	Contact customer service for repair.
Battery drains quickly	Battery failure	Contact customer service for repair.
Long charging time	Output current too low	Use an adapter with DC 5V 1A or higher.

TWO-YEAR WARRANTY

SEYSSO provides a two-year warranty from the date of purchase.

The warranty covers defects related to manufacturing faults or improper product assembly.

Keep your proof of purchase for the duration of the warranty.

You can file a warranty claim if the product has been used properly and as intended.

Warranty Does NOT Cover:

1. Surface damages from normal use.
2. Damages caused by:
 - a) Accidental drops, impacts, misuse, commercial use, neglect, or device modifications.
 - b) Improper maintenance or lack of care.
 - c) Connection to an incorrect power source.
 - d) Repairs or attempted repairs by unauthorized personnel.

Any attempt to open or dismantle the device, its parts, or accessories will void the warranty.

Warranty claims must be supported by proof of purchase and must be made within the warranty period.

Note: The handle material is delicate—handle with care. Scratches or wear on the surface are considered normal use damage and are not covered by the warranty.

The warranty covers functional components affecting device performance.

The instructions apply to model number: SE.WF.010.BLK.1, SE.WF.010.WHT.1

MANUFACTURER:

WM Brands sp. z o.o.
ul. Dziupli 12
02-454 Warsaw, Poland
Made in P.R.C.

CONTACT:

customer service: seysso@seysso.com
complaint service: reklamacje@seysso.com
product safety: kontakt@seysso.com
www.seysso.pl

Water resistance (handle): IPX7. Protection against immersion in water up to a depth of 1 meter for 30 minutes.

The irrigator has a built-in lithium-ion battery.

The product complies with the requirements of the European Union „New Approach“ directives.

The packaging is recyclable.

The instruction manual is recyclable.

This device is marked with the symbol of a crossed-out waste bin in accordance with the Polish Act on Waste Electrical and Electronic Equipment and the WEEE Directive 2012/19/EU. This marking informs that the device, after the end of its service life, must not be disposed of with other household waste. The user is obligated to return it to an electrical and electronic waste collection facility. Collection facilities, including local collection points, stores, and municipal units, create a system enabling the return of this device. Proper handling of electrical and electronic waste helps prevent harmful effects on human health and the environment caused by the presence of hazardous components and improper storage and processing of such equipment.

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